



eRad

Radiology Reporting System

Complete User Manual

Version 1.0.4 · March 2026

For radiologist, MRT, clinics, Hospitals and Diagnostic imaging centers

Table of Contents

1. Getting Started	4
1.1 System Requirements	4
1.2 First-Run Setup Wizard	4
Step 1 — Clinic Information	5
Step 2 — Admin Account.....	5
1.3 Logging In	6
2. Interface Overview	10
2.1 Navigation Sidebar	10
2.2 Dashboard	10
2.3 Theme Customization	11
3. Creating Reports	12
3.1 Patient Information	12
Searching for an Existing Patient.....	12
3.2 Selecting Exam & Service	12
3.3 Using the Rich Text Editor	13
3.4 Applying Templates.....	13
Saving a Template from Current Content	13
3.5 Saving & Finalizing	13
Mobile Floating Menu	14
4. Print Preview & Printing.....	15
4.1 Print Preview	15
4.2 Printing Reports	15
4.3 Exporting to PDF.....	15
5. Patient Records	16
5.1 Searching Records	16
5.2 Patient History	16
5.3 Record Actions	17
6. Report Templates.....	18
6.1 Browsing Templates	18
6.2 Creating a Template.....	18
7. Settings	20
7.1 Clinic Profile.....	20
7.2 Print Settings — Page Setup.....	20
7.3 Print Settings — Header	20
Per-Line Typography	21
7.4 Print Settings — Footer	21
7.5 Signature.....	22
8. User Management.....	24
8.1 User Roles.....	24
8.2 Adding Users	24

- 8.3 Digital Signatures 25
- 8.4 Changing Passwords 25
- 9. Services & Pricing 26
 - 9.1 Managing Services 26
 - 9.2 Categories 26
- 10. License & Activation 28
 - 10.1 License Plans 28
 - 10.2 Activating a License Key 29
 - 10.3 License Transfer to a New Computer 29
 - On the old PC 29
 - On the new PC 29
- 11. Backup & Restore 30
 - 11.1 Local Backup 30
 - What Is Backed Up? 30
 - 11.2 Restoring Data 30
 - 11.3 Recommended Backup Schedule 30
- 12. Troubleshooting & FAQ 31
 - Q: The app shows a blank white screen on startup. 31
 - Q: I forgot my admin password. 31
 - Q: "This license key is already activated on another computer." 31
 - Q: "Cannot connect to license server" during activation. 31
 - Q: Print preview is blank or shows only the header 31
 - Q: Lists (bullets/numbering) look correct in the editor but not in print. 31
 - Q: Images in the header (logo) are not showing in the preview. 31
 - Q: The app is slow when searching records. 31
 - Q: How do I add a second radiologist? 31
 - Q: Can two users log in at the same time? 31
 - Q: Where is my data stored? 32
 - Q: How do I update eRad to a new version? 32
- Getting Support 32


1. Getting Started

Installation, first-run setup, and logging in

1. Getting Started

1.1 System Requirements

Field / Option	Description
Operating System	Windows 10 / Windows 11 (64-bit)
RAM	4 GB minimum, 8 GB recommended
Disk Space	500 MB for installation + space for reports and backups
Display	1024 × 700 minimum resolution (1280 × 800 recommended)
Internet	Required only for license activation — fully offline after that
Printer	Optional — any Windows-compatible printer or PDF printer

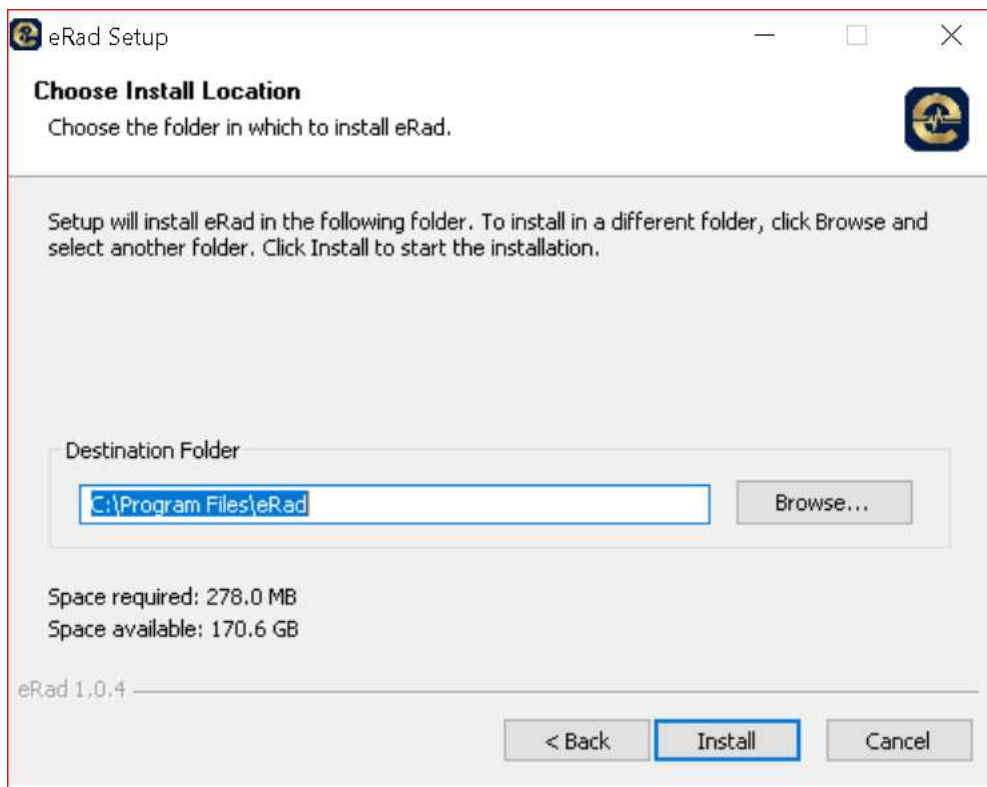

1.2 First-Run Setup Wizard

When eRad launches for the first time, a two-step setup wizard guides you through creating your clinic profile and the administrator account.

Step 1 — Clinic Information

Field / Option	Description
Clinic / Hospital Name *	Your clinic's primary English name. Appears on every printed report.
Address	Full street address of your clinic.
Phone	Main phone number (e.g. +251 9XX XXX XXXX).
Email	Contact email address.
Paper Size	Choose A4 (international) or Letter (US). Can be changed later in Settings.

Step 2 — Admin Account

Field / Option	Description
Full Name *	The administrator's full name (e.g. Dr. Solomon).
PRC License Number	Professional license number — optional, appears on printed reports.
Username *	Used to log in. Minimum 3 characters, lowercase.
Password *	Must be at least 6 characters. Case-sensitive.
Confirm Password *	Re-enter the password to confirm.

Note

Write down your admin username and password and store it safely. There is no "forgot password" option — if lost, data recovery may be difficult.

✓ Clinic Profile — **2** Admin Account — 3 Done

2 Create Admin Account


This is the main administrator account

Full Name *


PRC License Number

Username *

Password * Confirm Password *

 **Remember your password!**
If you forget it, data recovery may be difficult.


← Back ✓ Complete Setup



Welcome to eRad!


Let's set up your system in 2 quick steps

✓ Clinic Profile ——— ✓ Admin Account ——— 3 Done



Setup Complete!

X Hospital is ready.
Log in with your admin credentials to get started.

 [Go to Login →](#)


1.3 Logging In

After setup is complete you are directed to the Login screen. Enter your username and password, then press Enter or click Sign In.

Tip

Click the eye icon (👁) next to the Password field to show/hide your password as you type.

The login screen also shows your license status and expiry date at the bottom of the card.

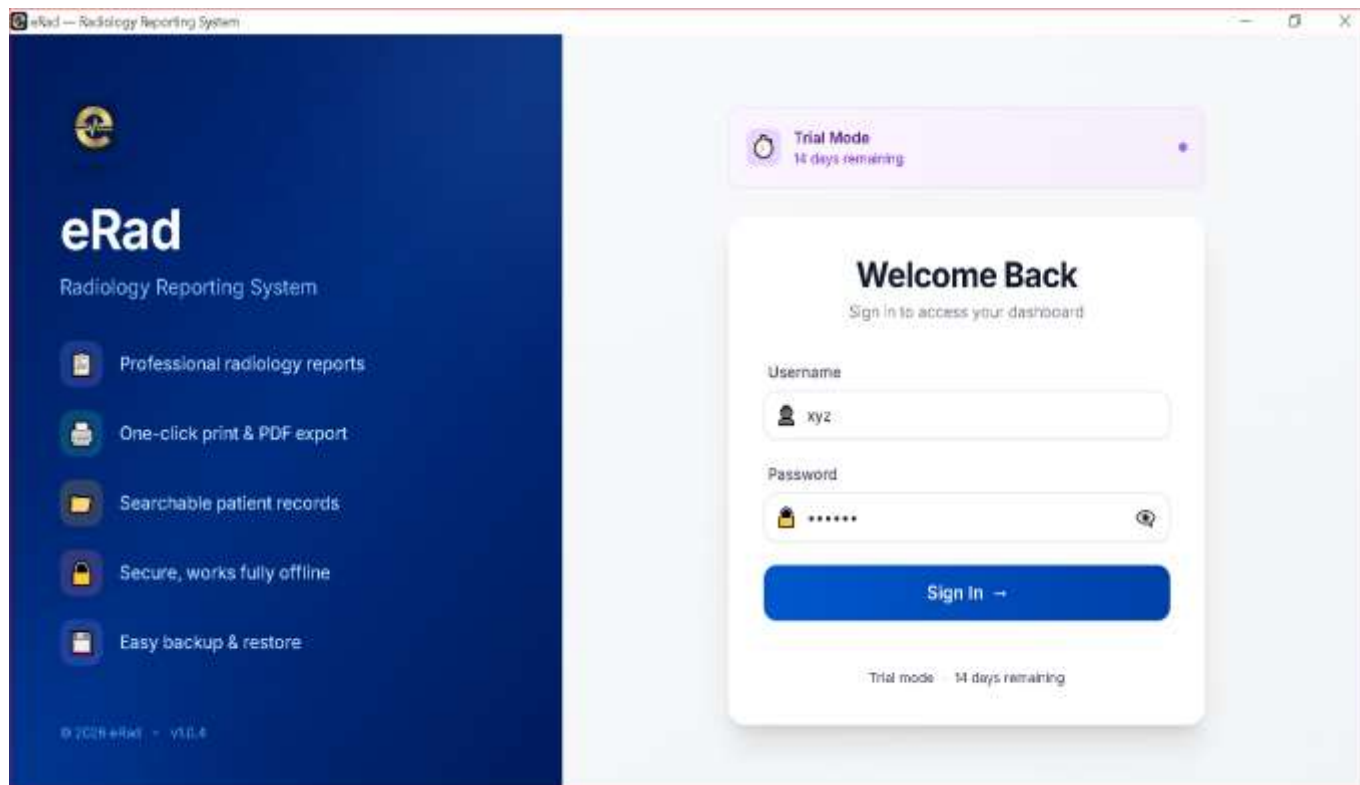

2. Interface Overview

Navigating the main application window

2. Interface Overview

2.1 Navigation Sidebar

The left sidebar is your main navigation hub. It shows your clinic logo, app version, and your role. Each menu item is role-restricted — you only see what you have access to.


Menu Item	Icon	Who Can Access
Dashboard	☐	All roles
New Report	☐	Admin, Radiologist, MRT
Records	☐	All roles
Templates	☐	Admin, Radiologist, MRT
Users	☐	Admin only
Services	☐	Admin only
Revenue	☐	Admin only
Settings	☐	Admin only

Tip


On smaller screens the sidebar collapses. Tap the blue circle button at the bottom-left of the screen to toggle it open.

2.2 Dashboard

The Dashboard is your home screen. It shows at-a-glance statistics including total reports, finalized reports, draft reports, patients seen today, and recent activity. Clicking any stat card navigates you to the relevant section.

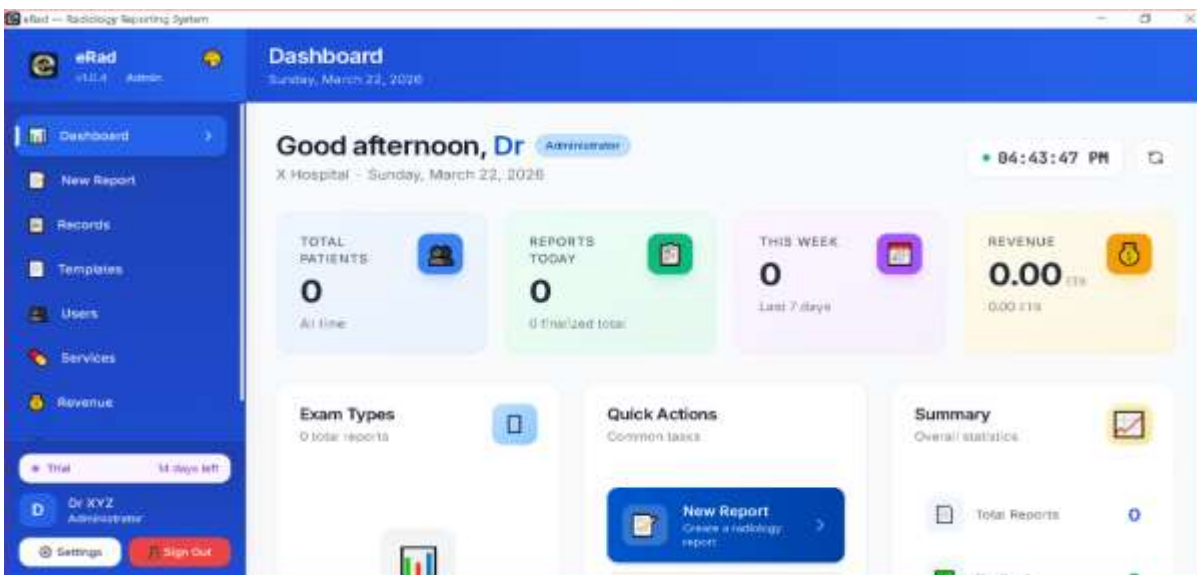

The "New Report" shortcut button is also available at the top-right of the Records page header for quick access.

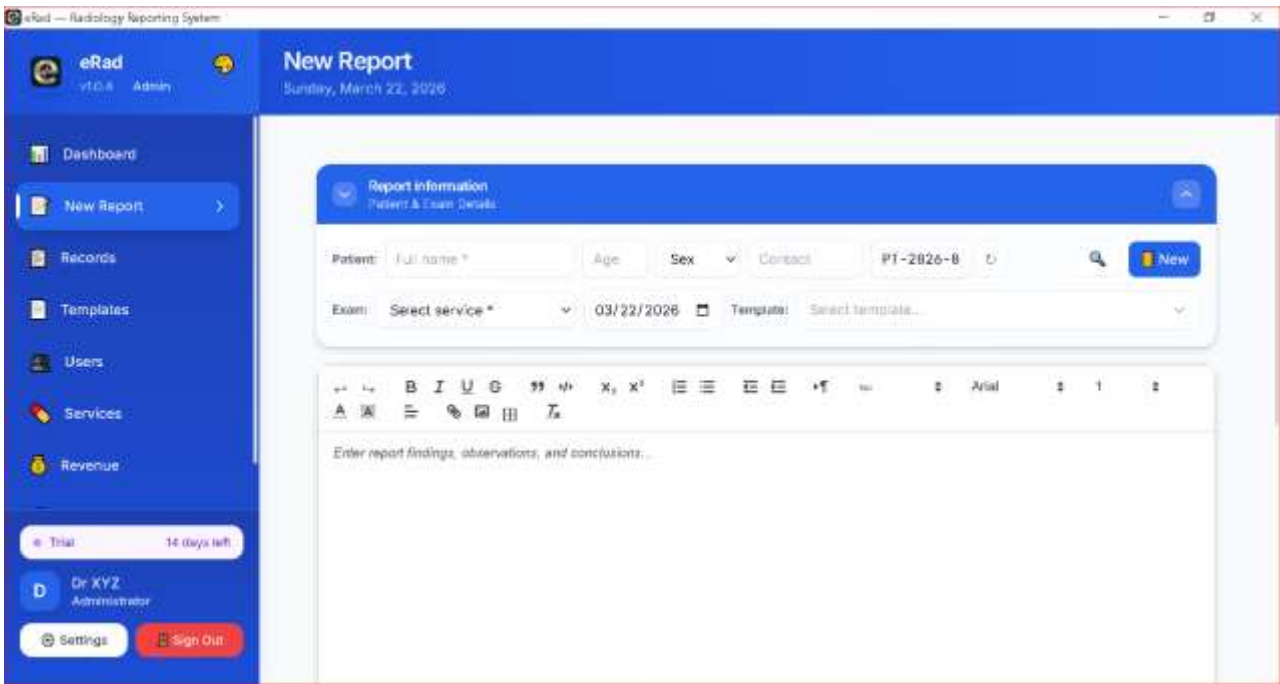

2.3 Theme Customization

eRad includes 5 color themes for the sidebar and header. Click the  paint palette icon in the top-right corner of the sidebar to open the theme picker.

Field / Option	Description
Professional Blue	Default — rich blue sidebar (recommended for clinical environments)
Corporate Slate	Neutral slate-gray — professional and subtle
Neutral Gray	Dark gray — low-contrast, easy on the eyes
Medical Green	Emerald green — classic medical aesthetic
Clean White	Light sidebar with dark text — high contrast

 **Tip**

Your chosen theme is saved automatically and remembered on next login.

3. Creating Reports

Step-by-step guide to writing and saving a radiology report

3. Creating Reports

Navigate to New Report from the sidebar. The page is split into two areas: the Report Information form at the top, and the Rich Text Editor below.

3.1 Patient Information

The patient section is inside the collapsible blue header bar. Click the bar to expand or collapse it at any time.

Field / Option	Description
Full Name *	Patient's full name. Required before saving.
Age *	Patient age in years. Numeric only.
Sex *	Male or Female. Required.
Contact	Phone number or other contact info. Optional.
Patient Code	Auto-generated unique code. Click ↻ to regenerate if needed.

Searching for an Existing Patient

Click the button to open the patient search panel. Type at least 2 characters to search by name, patient code, or phone number. Click a result to auto-fill all patient fields.

Tip

When an existing patient is selected, an "Existing" green badge appears. Click ✕ to clear and start a new patient entry.

3.2 Selecting Exam & Service

Field / Option	Description
Exam / Service *	Select from the grouped service dropdown. This sets the exam type and category used for billing and template filtering.
Date of Exam *	Defaults to today. Change if the exam was on a different date.
Template	Optional — pick a template to pre-fill the editor content. Search by name in the dropdown.

3.3 Using the Rich Text Editor

The editor supports full rich text formatting. The toolbar provides the following controls:

Field / Option	Description
Bold / Italic / Underline	Standard text formatting
Font family & size	Choose from a range of professional fonts and sizes
Text alignment	Left, center, right, or justify
Bullet lists	Unordered bullet points
Numbered lists	Ordered numbered lists
Tables (☐)	Insert a table — hover over the grid to choose rows and columns, or enter custom dimensions. Right-click a cell for row/column options.
Indent / Outdent	Increase or decrease list indentation
Undo / Redo	Standard undo/redo support

Tip

You can paste tables directly from Microsoft Word or Google Docs — eRad automatically cleans the formatting while preserving table structure.

3.4 Applying Templates

Templates pre-fill the editor with standard report text. After selecting a template, the content appears in the editor where you can edit it freely. Placeholder tokens are replaced automatically with patient data.

Placeholder	Replaced With
{patient_name}	Patient's full name
{age}	Patient's age
{sex}	Patient's sex (Male / Female)
{date}	Date of exam
{clinical_indication}	Clinical indication field
{technique}	Technique field

Saving a Template from Current Content

If you have written report content you want to reuse, click the purple Save button next to the template selector. Give the template a name and confirm. To update an existing template with your edits, click the amber Update button (only visible when a template is selected).

3.5 Saving & Finalizing


Field / Option	Description
Save Draft	Saves progress. The report remains editable. A draft badge appears.
Finalize	Marks the report as final. Appears with a green Finalized badge. Finalized reports can still be edited by authorized users.
Print	Sends the saved report to the printer. Requires the report to be saved first.

Field / Option	Description
Print Preview	Opens a visual preview of how the report will look when printed.

□ **Note**

Save your report as a Draft first before printing or finalizing. Unsaved content cannot be printed.

Mobile Floating Menu

On smaller screens, a blue floating  button appears at the bottom-right. Tap it to access Save Draft, Finalize, Print, Preview, and New Report actions.

4. Print Preview & Printing

Previewing, printing, and exporting reports to PDF

4. Print Preview & Printing

4.1 Print Preview

Click Print Preview (□) to open a full-page simulation of the printed report. The preview renders all pages accurately, showing exactly how the header, footer, and body content will appear on paper.

The preview dynamically measures header and footer heights so that body content never overlaps them — even on multi-page reports. Each page is separated by a gray "— Page N of M —" label for easy navigation.

□ **Tip**

You can print directly from the preview window using the blue □ Print button in the top-right corner of the modal.

4.2 Printing Reports

Click Print (□) in the report editor toolbar. The Windows print dialog opens where you can select your printer, paper size, and number of copies. The report must be saved at least once before printing.

4.3 Exporting to PDF

From the Records page, select a report and click PDF in the action panel on the right. A Save As dialog will appear — choose a location to save the PDF file.

□ **Tip**

PDF export uses the same layout engine as the print function, so the output is pixel-perfect with your print settings.

5. Patient Records

Searching, viewing, and managing patient reports

5. Patient Records

5.1 Searching Records

The Records page shows the 3 most recent reports by default. Use the search bar and filters to find specific records.

Field / Option	Description
Search bar	Search by patient name, report number, or exam type. Results update as you type (300 ms debounce).
Status filter	Filter by All, draft, or final.
Date From / To	Filter by exam date range.
View All Records	Click to load all records with full pagination (10 per page).

5.2 Patient History

Every patient row has a "History" link below their name. Click it to open the Patient History panel on the right sidebar, showing every visit for that patient.

From the Patient History panel you can click any past report to view its details, or click New Visit to create a new report for the same patient (their details are pre-filled automatically).

5.3 Record Actions

Click any report row to open its detail panel on the right. Available actions depend on your role:

Field / Option	Description
Edit Report	Opens the report in the New Report editor. Available to the report's author and admins.
View Content	Shows the raw report text in a modal without headers/footers.
Print	Sends the report to the printer.
PDF	Exports the report as a PDF file.
Finalize Report	Marks a draft report as final. Available to authorized users.
Delete Report	Permanently deletes the report. Admin only. Cannot be undone.

6. Report Templates

Creating and managing reusable report templates

6. Report Templates

Templates are pre-written report structures organized by exam category. They save time by pre-filling the report editor with standard content that you customize per patient.

6.1 Browsing Templates

The Templates page has a sidebar showing all exam categories. Click a category to filter the list. Use the search box to find templates by name. Click a template to preview its content on the right.

The navigator at the bottom of the sidebar lets you page through templates one at a time, or use the dropdown to jump directly to a specific template.

6.2 Creating a Template

Click the + New button in the sidebar header. Fill in:


Field / Option	Description
Template Name *	A descriptive name, e.g. "Normal Chest X-Ray" or "Brain MRI — No Acute Findings".
Category *	Must match an existing service category so templates appear when that exam is selected.
Template Content *	The full report text. Use the rich text editor with placeholders.

Tip

Categories come from the Services page. If a category is missing, add the service first, then create the template.

Note

Deleting a template that is linked to existing patient reports is blocked. You must update or delete those reports first.

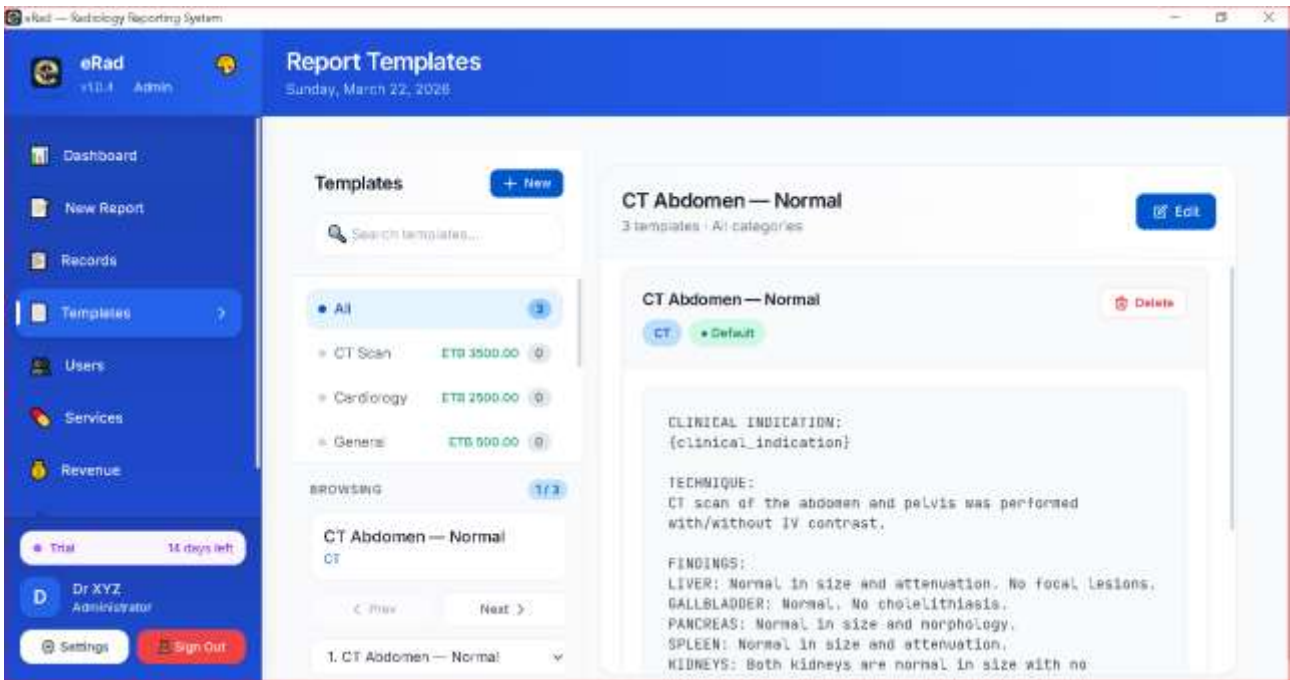

7. Settings

Clinic profile, print layout, and system preferences

7. Settings

Settings is accessible from the sidebar (☐) and is restricted to Admin users. It is organized into five tabs: Clinic Profile, License, Backup, Print, and System.

7.1 Clinic Profile

Configure your clinic's identity. This information appears on every printed report.

Field / Option	Description
Primary Name (English) *	Main clinic name, e.g. "St. Luke's Radiology Center".
Name in Local Language	Secondary name line in Amharic, Arabic, or any language. Displayed below the primary name.
Slogan / Tagline	Optional short tagline shown in the header.
Address, Phone, Fax, Email, Website	Contact details printed on every report.
Tax ID / VAT No.	Optional tax identification number.
Ministry / Registration No.	Optional regulatory registration number.
Default Paper Size	A4, Letter, or Legal.
Custom Print Header Text	Extra accreditation or certification text below the clinic name.
Clinic Logo	Upload a PNG or JPG logo (recommended 300×100 px). Enable/disable in Print settings.

☐ Done

Click Save Clinic Profile after making changes.

7.2 Print Settings — Page Setup

Field / Option	Description
Paper Size	A4, Letter (US), or Legal.
Body Font Family	Font used for the report body text (12 professional options).
Body Font Size	Slider from 7 pt to 20 pt.
Left / Right Margin	Side margins in millimeters.
Body Bottom Margin	Bottom margin in millimeters.

7.3 Print Settings — Header

The page header is pinned to the top of every page. Body content starts immediately below it — there is never any overlap.

Field / Option	Description
Show Header	Toggle the header on or off globally.
Padding Top / Bottom	Space inside the header above the clinic name and below the divider line.
Show Clinic Logo	Toggle the logo on/off. Upload the logo in Clinic Profile.
Logo Position	Left, Center, or Right.
Logo Height	Slider from 20 px to 120 px.
Header Font Family	Separate font from the body, applies to all header lines.
Header Text Alignment	Left, Center, Right, or Justify — applies to all header text.

Per-Line Typography

Each line in the header has its own font size and color control:

Field / Option	Description
① Clinic Name (Primary)	Size and color of the main English clinic name.
② Local Language Name	Size and color of the Amharic / Arabic name.
③ Slogan / Tagline	Size, color, and italic toggle for the tagline.
④ Contact Line	Size and color for address · phone · fax · email · website.
⑤ Registration / Ministry No.	Size and color for the registration number.
⑥ Extra Header Lines	Size, color, and additional free-text line.

A Live Preview panel below the controls shows exactly how the header will look before you save.

7.4 Print Settings — Footer



The footer is pinned to the bottom of every page. Body content never overlaps it.


Field / Option	Description
Show Footer	Toggle footer visibility.
Padding Top / Bottom	Internal spacing above and below the footer text.
Footer Text	Main disclaimer or notice text shown on the left.
Footer Font Family	Font for the footer text.
Footer Font Size	Size in pt.
Footer Text Color	Color picker with preset swatches.
Footer Text Alignment	Left, Center, Right, or Justify.
Right Stamp	"Printed: [date] · eRad" shown on the right side. Has its own size and color controls.

A Footer Preview panel shows a real-time preview of the footer layout.

7.5 Signature

Enable or disable the radiologist signature block at the bottom of each report. Upload individual signatures from the Users page (Users → select user → Upload Signature).

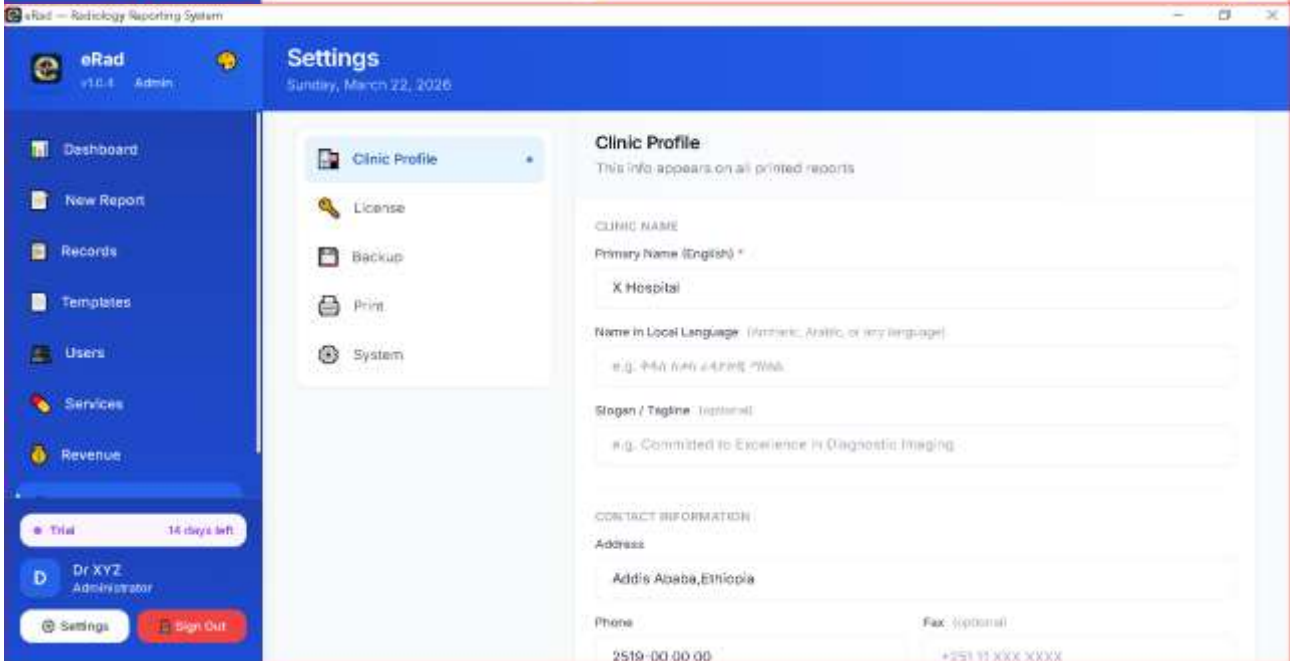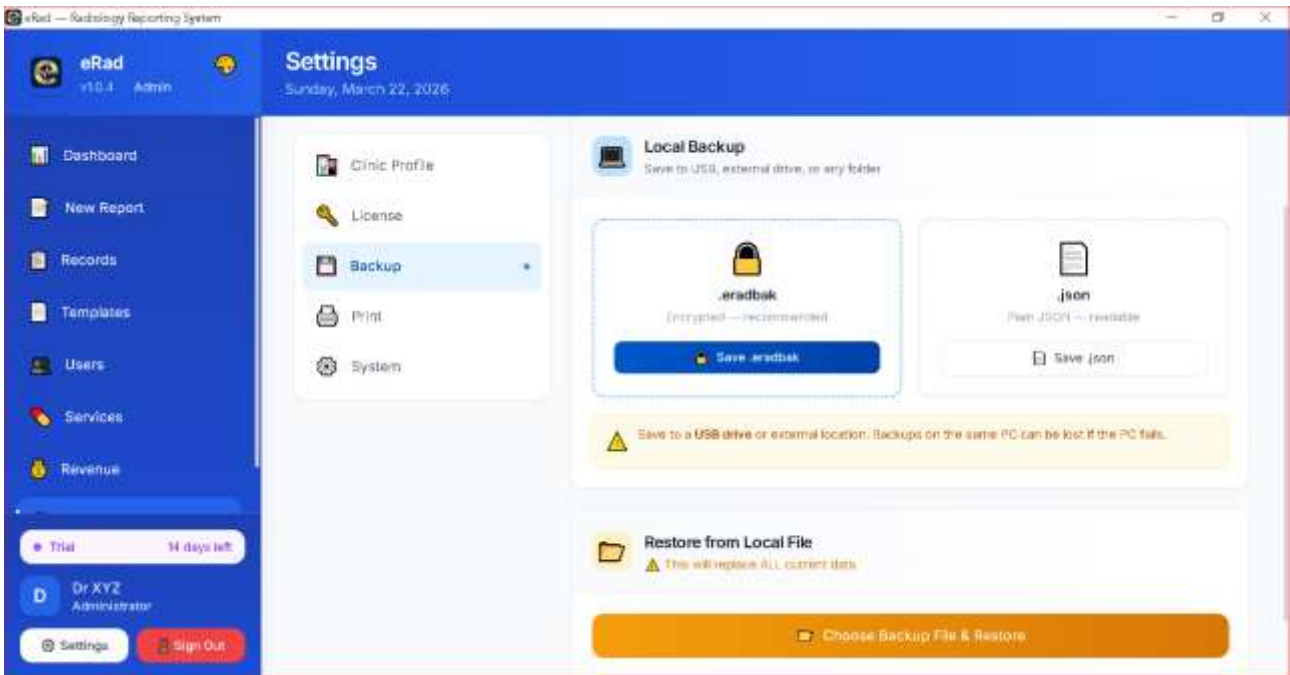

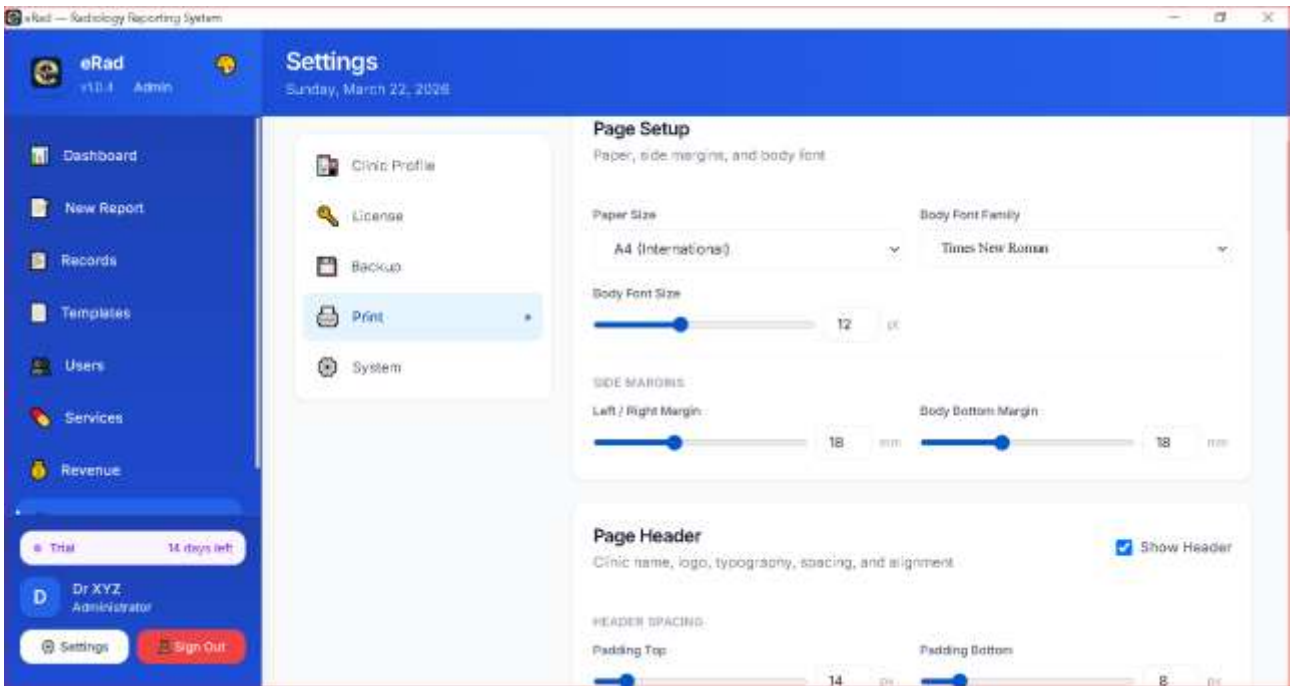

8. User Management

Adding users, assigning roles, and managing credentials

8. User Management

8.1 User Roles

eRad has four roles, each with different access levels:

Role	Label	Access Level
admin	<input type="checkbox"/> Admin	Full access — manage users, settings, backup, all reports
radiologist	<input type="checkbox"/> Radiologist	Create & finalize reports, manage templates
mrt	<input type="checkbox"/> MRT	Create & edit reports, register patients
receptionist	<input type="checkbox"/> Receptionist	Register patients only, view records

8.2 Adding Users

Go to Users in the sidebar. Click + Add User — a panel slides in on the right.

Field / Option	Description
Full Name *	User's display name.
Username *	Login username — unique, cannot be changed after creation.
Role *	Admin, Radiologist, MRT, or Receptionist.
PRC License No.	Professional license number. Appears on printed reports.
Print Title	Title shown under the user's name on reports, e.g. "Medical Radiologic Technologist" or

Field / Option	Description
	"Radiologist, MD". Auto-populated based on role.
Password *	Minimum 6 characters (new users only).

Tip
 To disable a user without deleting them, click the Deactivate button on their card. Inactive users cannot log in but their historical reports remain intact. You can reactivate them at any time.

8.3 Digital Signatures


Radiologists can have a signature image attached to their profile. It appears below the report content when printing, above the name and title line.

To upload a signature: go to Users, select the user, then click Upload Signature. Any image format is accepted. Recommended: a transparent-background PNG of the actual handwritten signature, 400x150 px.

Tip
 Enable or disable signature printing globally in Settings → Print → Signature.

8.4 Changing Passwords

Select a user in the Users page, then click Change Password. If you are changing your own password, you must provide your current password first. Admins can change any user's password without knowing the old one.

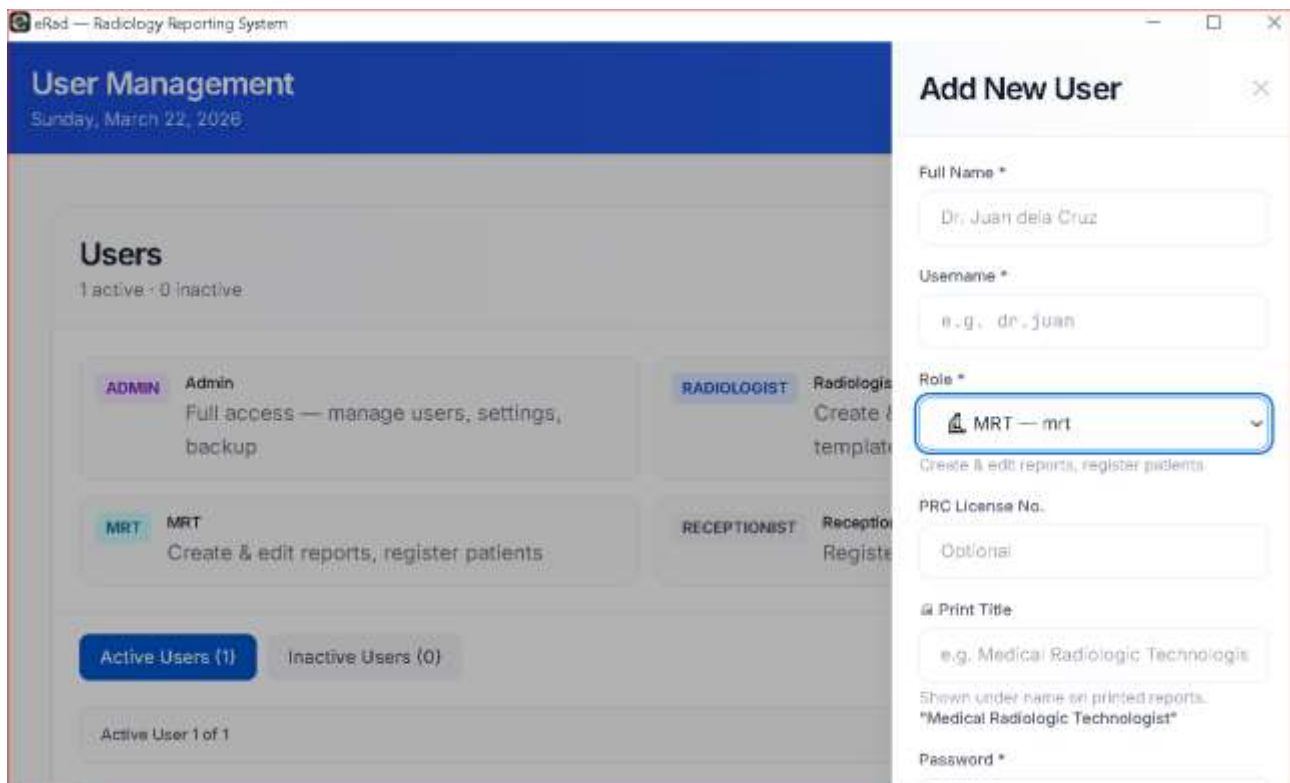

9. Services & Pricing

Managing your exam catalog and billing rates

9. Services & Pricing

Services are the exam types available in eRad. They are organized by category (X-Ray, MRI, CT, Ultrasound, etc.) and have a set price used for billing. eRad comes with a comprehensive default catalog of 90+ services.

9.1 Managing Services

Navigate to Services from the sidebar. Use the search bar and category filter to find services. Click + Add Service to create a new one.

Field / Option	Description
Service Name *	Descriptive name, e.g. "Brain MRI with Contrast".
Category *	X-Ray, MRI, CT, Ultrasound, ECHO, ECG, EEG, or Other.
Currency	ETB (Ethiopian Birr) or USD.
Price *	Numeric price used for billing.
Sort Order	Lower numbers appear first in dropdowns. Default 0.
Description	Optional short description.

Tip

Deactivating a service hides it from the New Report dropdown but preserves all historical billing records. Use "Show inactive" to view and restore them.

9.2 Categories


The service category is the link between services and templates. When a service is selected in a report, only templates with the matching category are shown in the template dropdown.

Tip

To add a new category, simply create a service with a new category name. The category will appear automatically throughout the app.

The screenshot shows the 'Service Catalog' page in the eRad system. The left sidebar contains navigation options: Dashboard, New Report, Records, Templates, Users, Services (selected), and Revenue. The main content area is titled 'Services & Pricing' and includes a '+ Add Service' button. Summary cards show 98 Active Services, 7 Categories, and an Average Price of Br 3586.22. A search bar and filter tabs (All, X-Ray, MRI, CT, Ultrasound, ECHO, ECG, EEG, Other) are present. The current view shows 'Page 1 of 98' with '98 services found'. A dropdown menu is open for '1. Brain CT (Non-contrast) (CT)'. The service details for 'Brain CT (Non-contrast)' are displayed, showing a price of Br 3500.00 and an 'Active' status.

The screenshot shows the 'Revenue & Billing' page in the eRad system. The left sidebar is identical to the previous page, with 'Revenue' selected. The main content area is titled 'Revenue & Billing' and includes a 'Refresh' button. A 'Quick Range' section allows filtering by 'Today', 'This Week', 'This Month', and 'This Year'. The current filters are 'This Month', 'From 02/28/2026', 'To 03/22/2026', and 'Radiologist: All Staff'. Four summary cards are displayed: 'TOTAL REVENUE 0.00 Br' (0 billed reports), 'COLLECTED 0.00 Br' (0% collection rate), 'OUTSTANDING 0.00 Br' (Awaiting payment), and 'TOTAL DISCOUNTS 0.00 Br' (Given across all reports). At the bottom, there are tabs for 'Overview', 'By Service', 'By Radiologist', and 'Recent Transactions'.



XYZ HOSPITAL

ኤክስ ዋይ ኤድ ሆስፒታል

Addis Ababa, Ethiopia

RD-2026-00021

2026-03-22

Patient Name	Mr X	Patient Code	PT-2026-0037
Age / Sex	55 yrs / Male	Date of Exam	2026-03-22

BILATERAL CAROTID DOPPLER ULTRASOUND

AOD(cm)	<4cm	N	MPA(cm)	0.9-2.9	N
LAD(cm)	3cm-4cm	N	MVA	4-6cm ²	N
LVD(cm)	4.3cm-5.9cm	N	EF	>55	60
IVS(cm)	2.6cm-4cm	N	E/A	0.8-2	N
IVSD(cm)	0.6cm-1.3cm	N	E DT(ms)	150-250	N
PWD(cm)	0.6cm-1.2cm	N	LVOT Vp(m/s)	<1.3	N
RA(cm)	2.9cmx4.5cm	N	AO Vp(m/s)	0.9-1.7	N
RVD(cm)	<4.2cm	N	TR Vp(m/s)	1.7-2.3	N
RVOT(cm)	1.8cm-3.8cm	N	PA Vp(m/s)	0.5-1.0	N
IVC	<1.8cm,>50%		E'(m/s)	0.8-1.3	N

SUMMARY

The left ventricle has normal cavity dimensions, wall thickness contractility, it has normal global systolic function and EF is 60 % and FS 33%; there is no sign of regional wall motion abnormality

- The right ventricle has normal cavity dimension and systolic function with TAPSE of 20mm
- The right and left atrium have normal dimension. There is no sign of tricuspid regurgitation.
- The mitral valve leaflets have normal morphology, mobility and valve opening.
- Normal E and A waves demonstrated on TDI.
- The aortic and pulmonic valves have normal leaflet morphology, mobility and valve opening; they have normal flow pattern on Doppler evaluation; there is no sign of stenosis or regurgitation.

CONCLUSION:

Normal Echocardiography

Dr XYZ
 MD, Radiologist

10. License & Activation

Understanding license plans and activating your software

10. License & Activation

10.1 License Plans

Plan	Code	Duration	Typical Price (ETB)
Trial	T	14 days	Free
Monthly	M	30 days	250
Quarterly	Q	3 months	600
Biannual	H	6 months	1,000
Annual	A	1 year	1,900
Lifetime	L	Forever	5,000

Tip

Paid license keys stack — activating a new key before expiry extends your current license by the new key's full duration. For

example, activating an Annual key 30 days before expiry gives you 365 + 30 days total.

10.2 Activating a License Key

Go to Settings → License → scroll to "Activate License". Type or paste your license key (format: ERAD-X-XXXXXXXX-XXXX) and click Activate.

Note

Internet access is required for activation. After successful activation, eRad works completely offline — no further internet connection is needed.

Contact your supplier via WhatsApp (+251 936 239 911) or email (davesolabay@gmail.com) to obtain a license key. Payment options include CBE, BOA, and Telebirr.

10.3 License Transfer to a New Computer

Each license is tied to one machine. To move eRad to a new PC, follow this checklist:

On the old PC

1. Go to Settings → License → Transfer to Another Computer.
2. Click Export License → save the .eradkey file to a USB drive.
3. Click Export All Data (Backup) → save the .eradbak file to the same USB drive.

On the new PC

4. Install eRad.
5. Go to Settings → License → Transfer → click Import License.
6. Select the .eradkey file from the USB drive.
7. Go to Settings → Backup → Restore → select the .eradbak file.

Done

After importing the license, eRad will restart. Then restore your backup. All patients, reports, and templates will be exactly as before.

11. Backup & Restore

Protecting your data with regular backups

11. Backup & Restore

eRad stores all data locally on your PC. Regular backups protect against accidental data loss, hardware failure, or Windows issues.

11.1 Local Backup

Go to Settings → Backup. Two backup formats are available:

Field / Option	Description
.eradbak (Encrypted)	Recommended. Encrypted binary format — safe to send or store externally.
.json (Plain)	Human-readable JSON. Useful for data inspection or migration.

Note

Save backups to a USB drive, external hard disk, or cloud storage (Google Drive, Dropbox). Backups saved only on the same PC are lost if the PC fails.

What Is Backed Up?

A backup includes ALL data: patients, reports, report templates, users, services, billing records, and settings. Nothing is excluded.

11.2 Restoring Data

Go to Settings → Backup → Restore. Click "Choose Backup File & Restore" and select a .eradbak or .json backup file.

Warning

Restoring PERMANENTLY overwrites all current data with the contents of the backup. Make a new backup of your current data before restoring if you might need it.

11.3 Recommended Backup Schedule

Field / Option	Description
Daily	High-volume clinics (20+ reports/day) — backup at end of each day.
Weekly	Moderate volume — backup at end of each working week.
Before upgrades	Always backup before installing a new version of eRad.
Before restore	Always backup current data before restoring from an old backup.

12. Troubleshooting & FAQ

Common questions and how to resolve them

12. Troubleshooting & FAQ

Q: The app shows a blank white screen on startup.

Wait 5–10 seconds for the loading spinner. If it persists, close and reopen eRad. If still blank, check that your Windows version is up to date.

Q: I forgot my admin password.

There is no self-service password reset. Contact your eRad supplier and provide your Machine ID (Settings → License). A database reset may be required.

Q: "This license key is already activated on another computer."

Each license key works on one machine only. Use the License Transfer workflow (Settings → License → Transfer) to move your license to the new PC. If the old PC is no longer available, contact your supplier.

Q: "Cannot connect to license server" during activation.

Check your internet connection. Disable any firewall or antivirus that may be blocking outbound HTTPS connections to erad.doxhall.com.

Q: Print preview is blank or shows only the header.

Save the report first (Draft or Finalize), then open Preview again. Unsaved content cannot be previewed.

Q: Lists (bullets/numbering) look correct in the editor but not in print.

This is a known quirk of some older system fonts. Try changing the Body Font Family in Settings → Print → Page Setup to Times New Roman or Georgia.

Q: Images in the header (logo) are not showing in the preview.

Make sure the logo file still exists at the path it was uploaded from. If the file was moved or deleted, re-upload it from Settings → Clinic Profile → Upload Logo.

Q: The app is slow when searching records.

eRad is optimized for fast search. If you have thousands of records, make sure your PC meets the minimum RAM requirement (4 GB). Close other heavy applications.

Q: How do I add a second radiologist?

Go to Users → Add User. Set the role to Radiologist. Each user has their own login, reports, and optional signature image.

Q: Can two users log in at the same time?

eRad is a single-user desktop application — only one user can be logged in per PC instance. Multiple PCs in the same clinic each need their own license.

Q: Where is my data stored?

All data is stored in C:\ProgramData\eRad\ on Windows. You can open this folder from Settings → System → Open Data Folder.

Q: How do I update eRad to a new version?

Download the new installer from erad.doxhall.com and run it. Your data is preserved. Always backup first before updating.

Getting Support

Field / Option	Description
WhatsApp	+251 936 239 911
Email	davesolabay@gmail.com
Website	https://erad.doxhall.com
Machine ID	Settings → License → Machine ID (always include when contacting support)

eRad Radiology Reporting System — User Manual

© 2026 eRad Radiology Systems. Document generated 22 March 2026.